

Unheard-of innovation.

# THE NEW EARBREEZE® DRY

We're all familiar with it: *the feeling of wet ears after showering, swimming or water sports.* Unfortunately, in many cases, this unpleasant sensation of pressure is often followed by inflammation in the ear canal, which can become chronic.

This is due to a hollow near to the eardrum, which gathers residual moisture and fosters the breeding of bacteria and germs. Hopping up and down on one leg won't help here, and pushing a cotton swab deep into your ear actually invites additional risks.

This is why we have developed earbreeze® dry: *the first practical aid in the world, an ear dryer that dries your ears easily and gently.*



earbreeze® dry  
Warm air ear dryer

## LISTEN TO YOUR EARS



Splash-proof  
ear case



Replaceable  
battery



Rechargeable  
200x

"Many of my patients suffer from recurring infections due to their wet ear canals at high levels of humidity, after water sports or during sweat-inducing activities. earbreeze® dry is able to prevent this and is safe and easy to use."

DR. PETER RAMBERGER,  
EAR, NOSE AND THROAT SPECIALIST,  
VIENNA

CURIOUS?




Find out more now.

earbreeze GmbH  
Wipplingerstraße 1/1.DG, 1010 Vienne | Austria  
earbreeze.com

earbreeze® dry  
Warm air ear dryer

## NO MORE WET EARS

showering • swimming • water sports

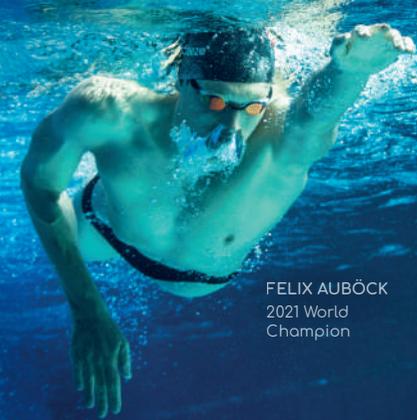

FELIX AUBÖCK  
2021 World  
Champion



The new earbreeze® dry gets rid of residual moisture in your ear canal swiftly and gently, providing effective prevention against ear infections.

earbreeze.com

## READY FOR EARBREEZE® DRY?

earbreeze® dry  
helps prevent:

any unpleasant feelings  
of pressure in the ear

earache and infections

chronic conditions such  
as swimmer's ear, diver's  
ear and surfer's ear

any build-up of moisture  
when wearing hearing  
aids



Developed in Austria.

INTER-  
NATIONALLY  
PATENTED

Patented "tongue"  
of the earbreeze® dry




The new earbreeze® dry was developed by experts, based on the latest scientific findings, to help those affected in their day-to-day lives. This is how the first gentle ear dryer with an internationally patented nozzle design was developed.

The core element of the nozzle design is a "tongue" for swirling the airflow at body temperature. This enables the airflow to gently reach the outer ear canal, and the air carrying the residual moisture flows through the earpiece to the outside in a controlled manner, protecting the eardrum.

## EASY, EASIER, EARBREEZE® DRY

**How the new earbreeze® dry works:**

You activate your earbreeze® dry with the push of a button. After a brief warming-up period to 37°C, you hear and see that your earbreeze® dry is ready.

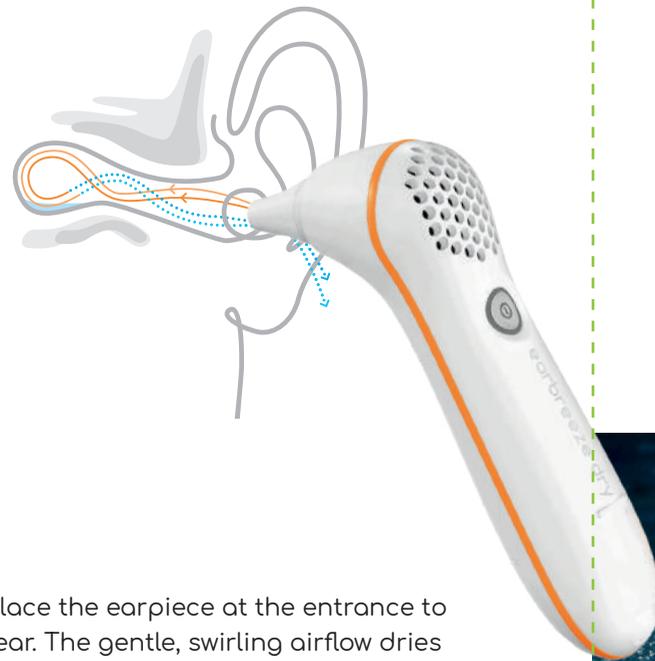

You place the earpiece at the entrance to your ear. The gentle, swirling airflow dries your outer ear canal at body temperature. One application takes 45 seconds and you can repeat it as many times as you need.

The right shape, every time.

## FOR BIG AND SMALL EARS ALIKE

There are two different sizes of earpiece to meet the needs of varying ear sizes and shapes. These are easy to replace and clean.

There are also three speed levels for the air generated by the device to suit your sensitivity. This makes the new earbreeze® dry helpful for the whole family.

Handy  
size



Integrated  
battery



Easy to  
transport



Sound all around.

## USE WITH HEARING AIDS

Moisture can also frequently build up with the use of hearing aids and in-ear headphones, due to the reduced ventilation of the ear canal. In both cases, you can use the new earbreeze® dry to increase your wearing comfort and prevent infections.

Simply briefly remove your hearing aid from your ear a few times a day and dry your ear canal using the earbreeze® dry. Then reinsert your hearing aid as you normally would.

